



## Run Info

Host Name	GXB02097 (localhost)
Position	X2
Experiment Name	VX06_H61211b
Sample ID	H61211b
Run ID	acfd916c-b527-4e2a-b6ce-049a28c60626
Acquisition ID(s)	db62ecd5accb5327d8d35601c7597e3a72426a1b, 3433ed21acc27b331e98aacba6f60fed133940c8
Flow Cell Id	FAL83665
Start Time	July 22, 10:53
Run Length	3d 0h 4m

## Run Summary

Reads Generated	10.72 M
Passed Bases	4.12 Gb
Failed Bases	2.78 Gb
Estimated Bases	7 Gb


## Run Parameters

Flow Cell Type	FLO-MIN106
Kit	SQK-DCS109
Initial bias voltage	-180 mV
FAST5 output	Enabled
FASTQ output	Enabled
BAM output	Disabled
Bulk file output	Disabled
Active channel selection	Enabled
Basecalling	Enabled
Specified run length	72 hours
FAST5 reads per file	4000
FAST5 output options	vbz_compress,fastq,raw
FASTQ reads per file	4000
FASTQ output options	compress
Mux scan period	1 hour 30 minutes
Reserved pores	0 %
Basecall model	High-accuracy basecalling
Read filtering	min_qscore=9


## Versions

MinKNOW	21.05.12
MinKNOW Core	4.3.7
Bream	6.2.5
Guppy	5.0.12

Cumulative Output Reads

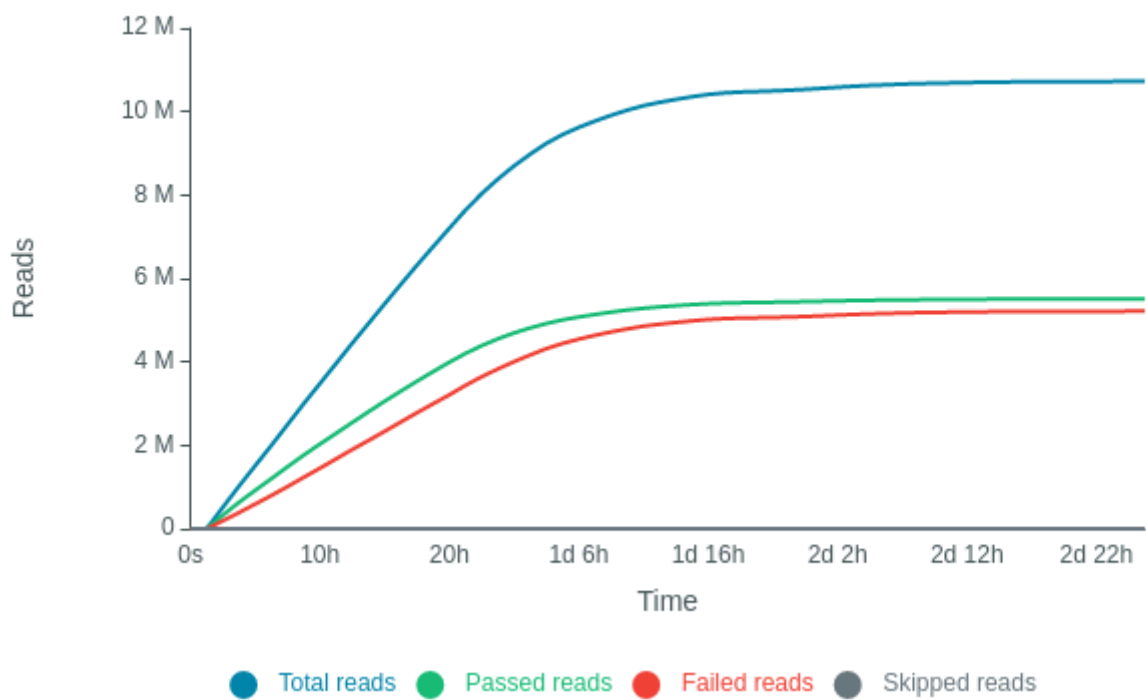


Cumulative Output Bases

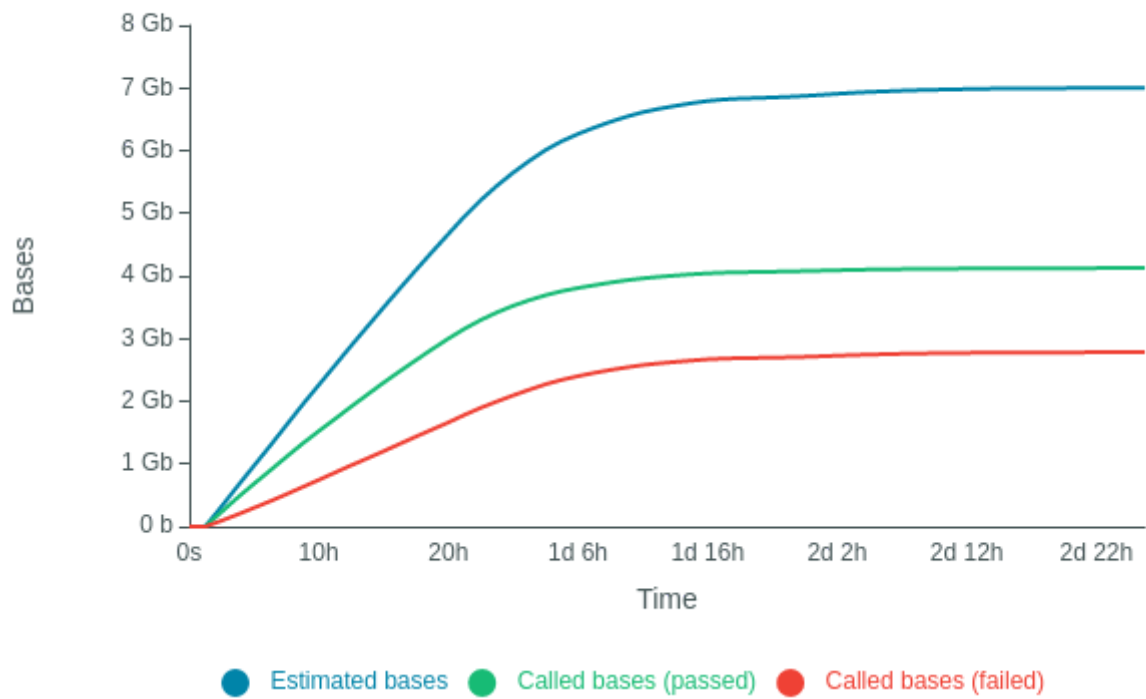


**Read Length Histogram Estimated Bases - Outliers Discarded**

Estimated N50: 785 b

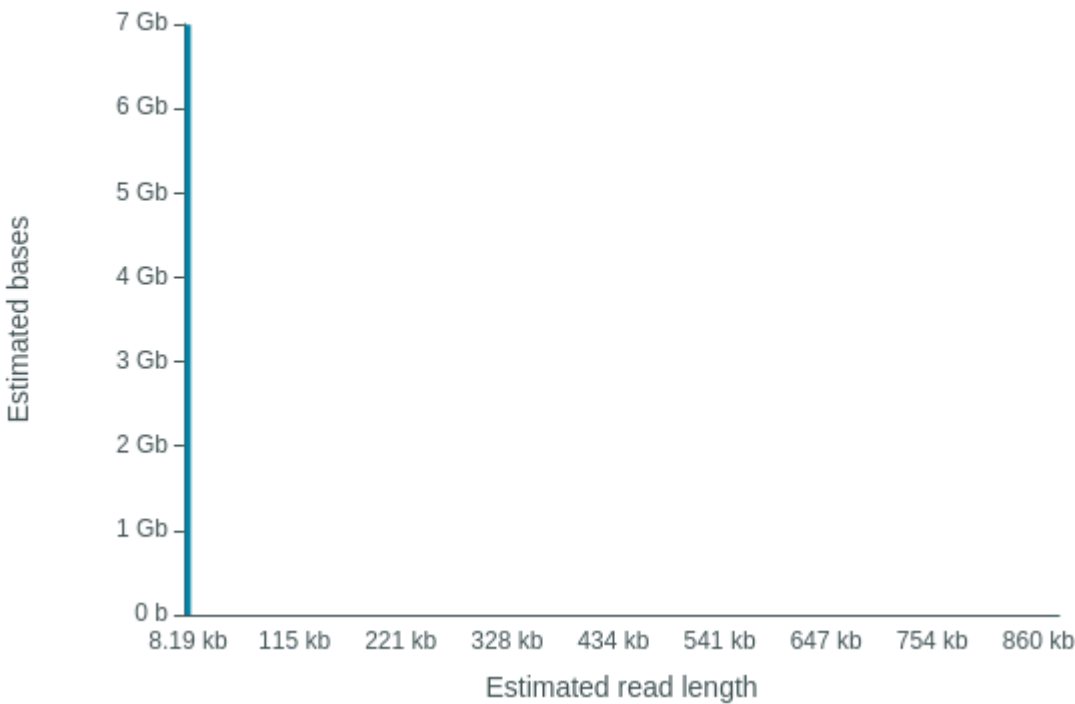


**Read Length Histogram Basecalled Bases - Outliers Discarded**

Estimated N50: 812 b

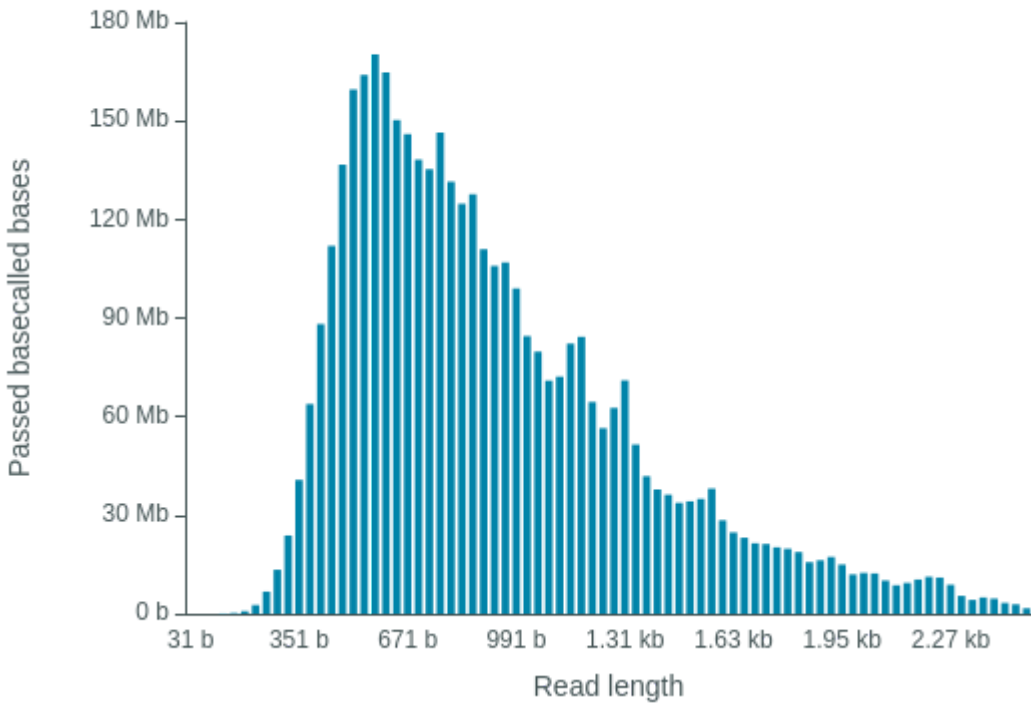


**Read Length Histogram Estimated Bases**

Estimated N50: 785 b

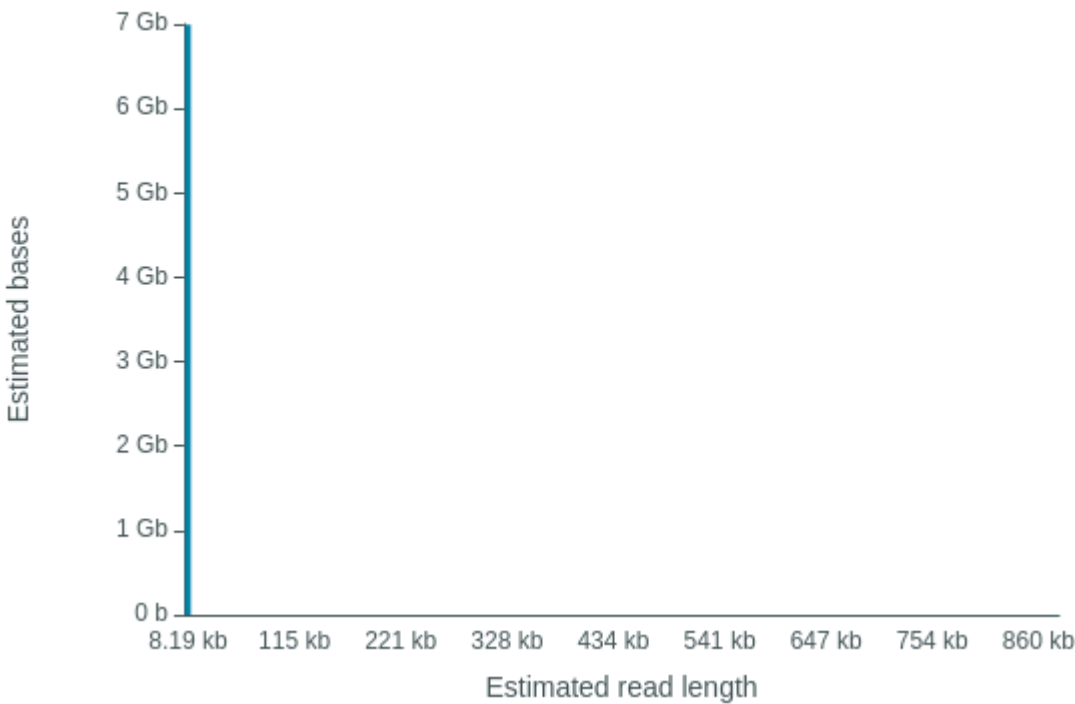


**Read Length Histogram Basecalled Bases**


Estimated N50: 812 b

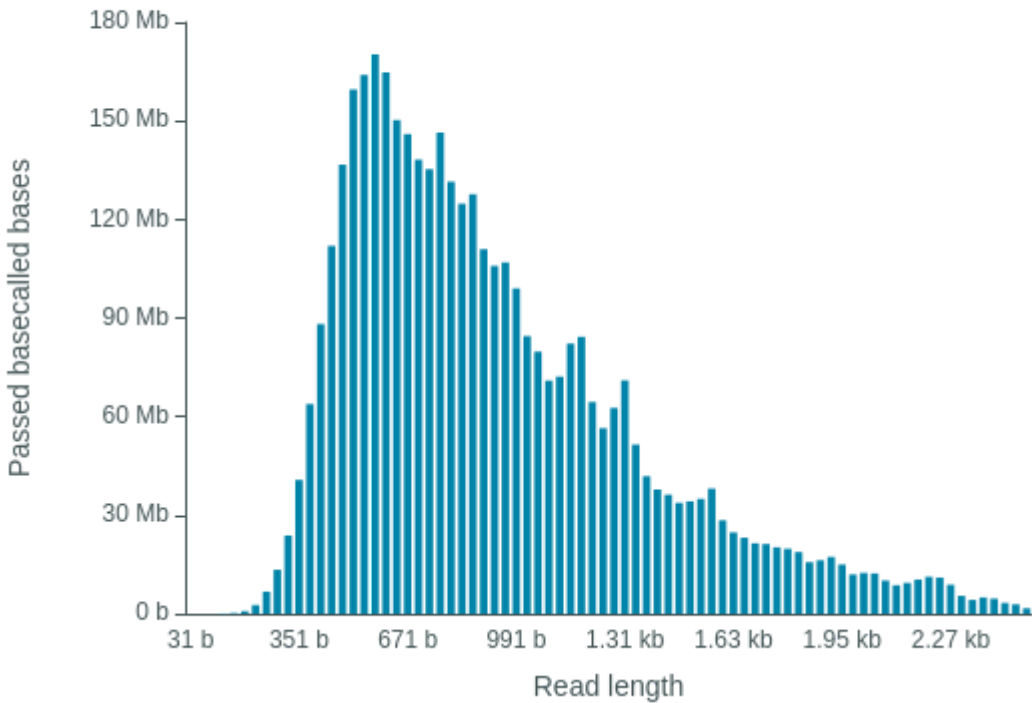


Duty Time Grouped

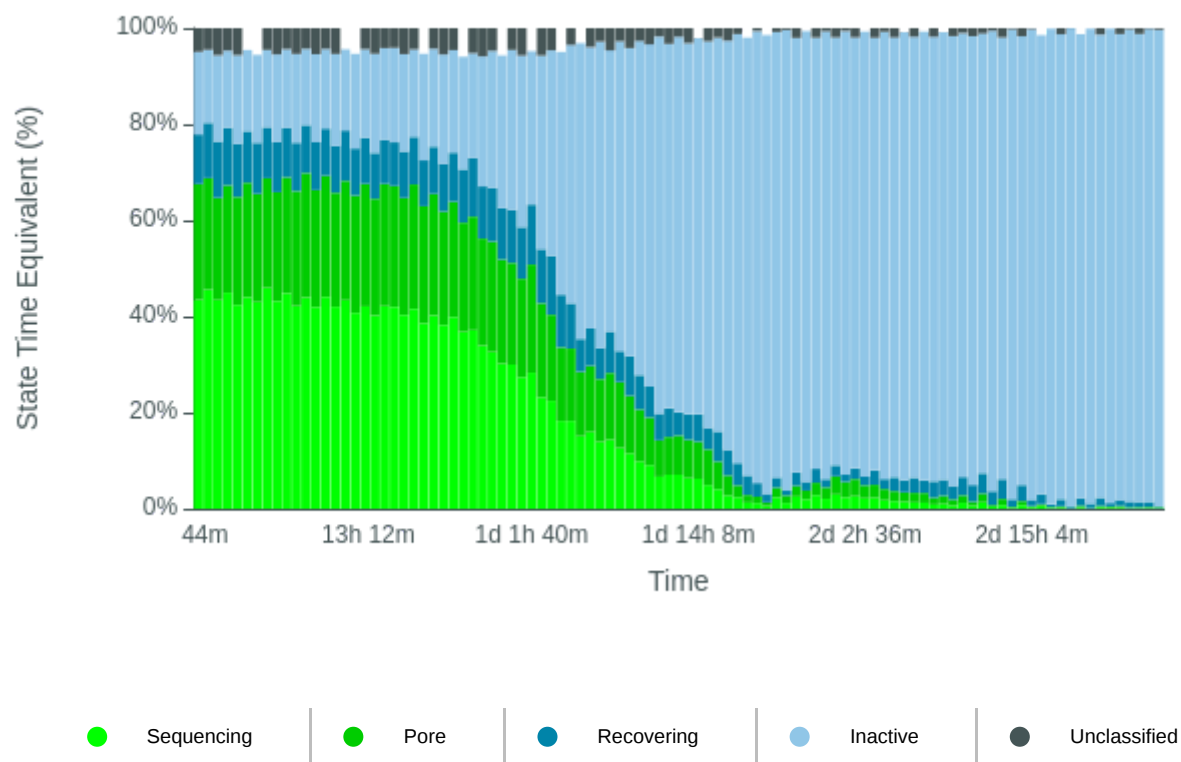


Duty time Categorised

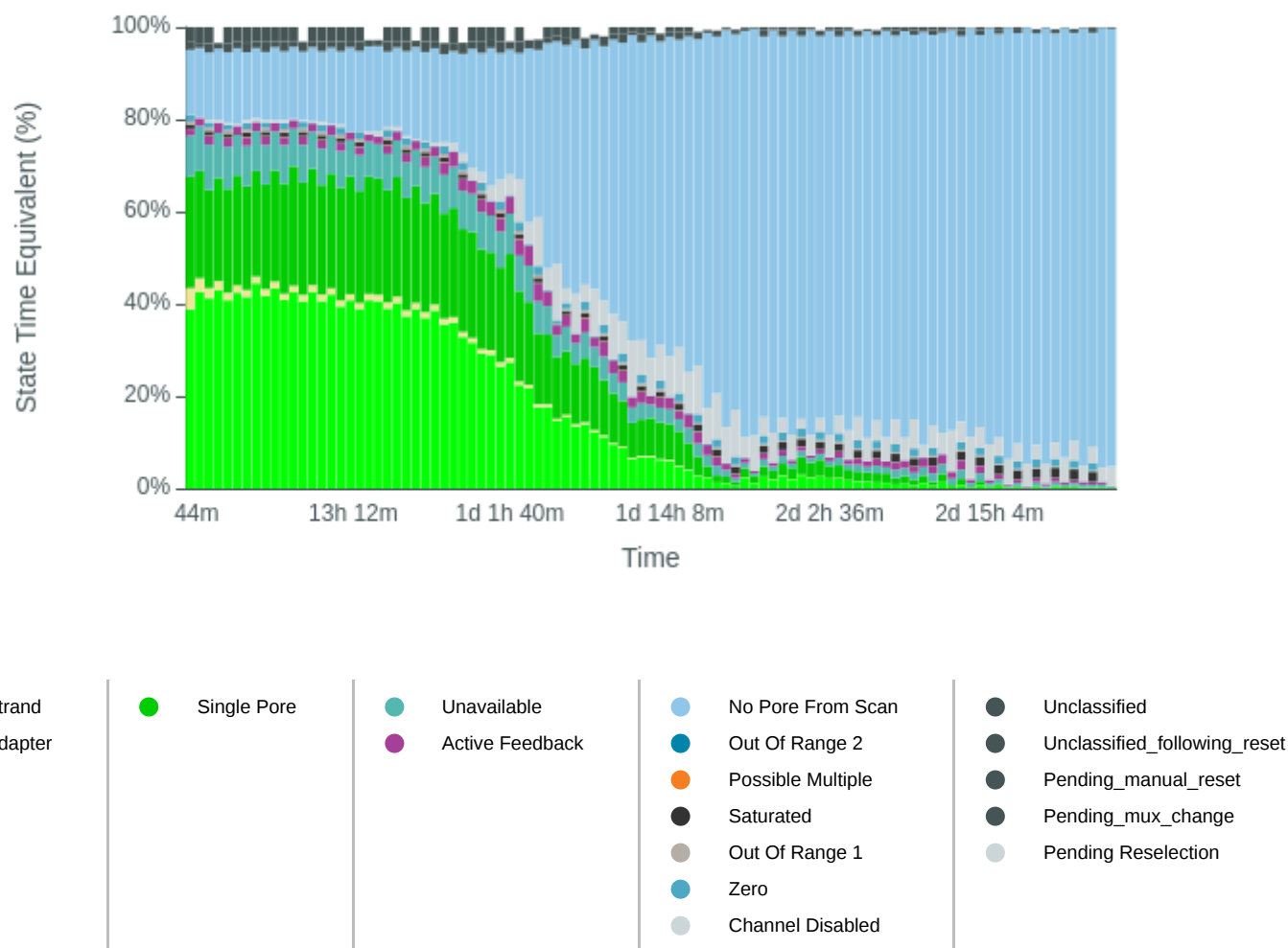


Mux Scan Grouped

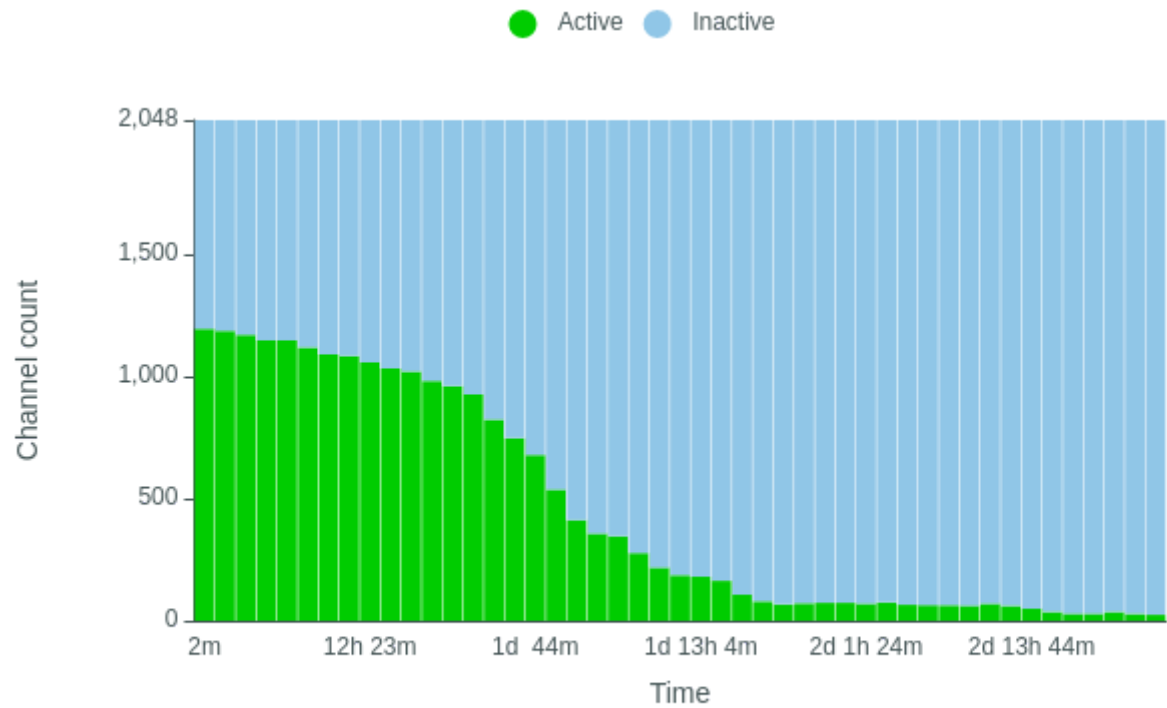


Mux Scan Categorised

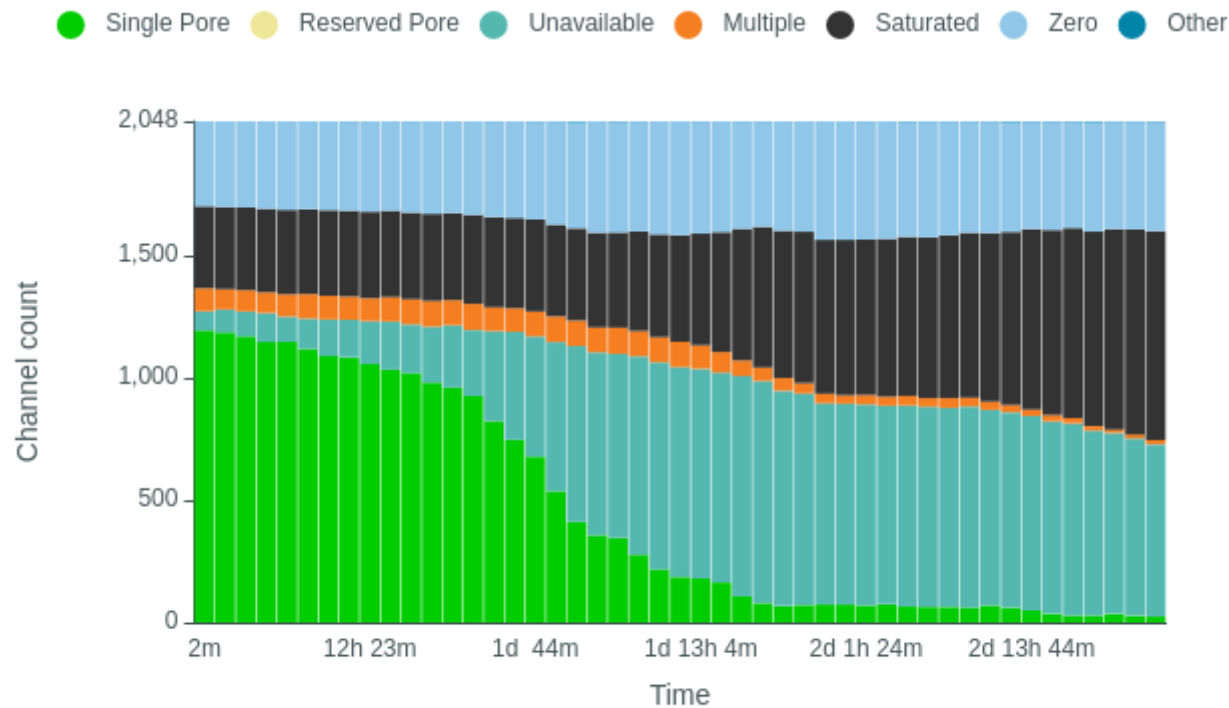


Temperature History

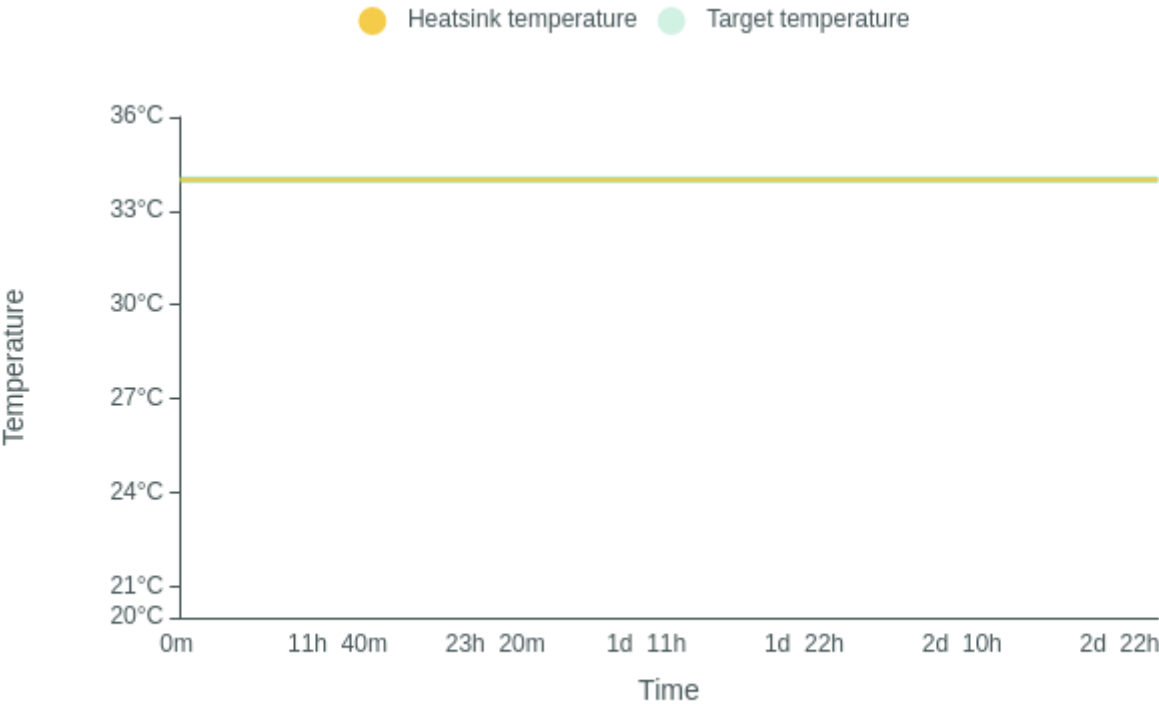


Bias Voltage History

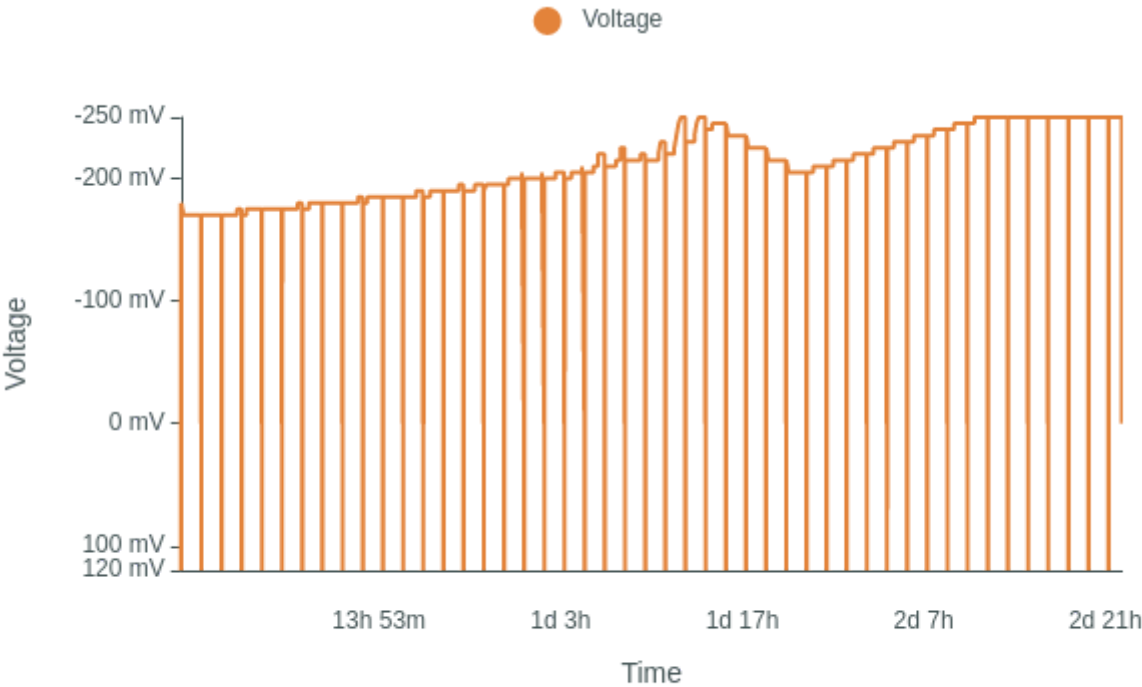

Translocation Speed

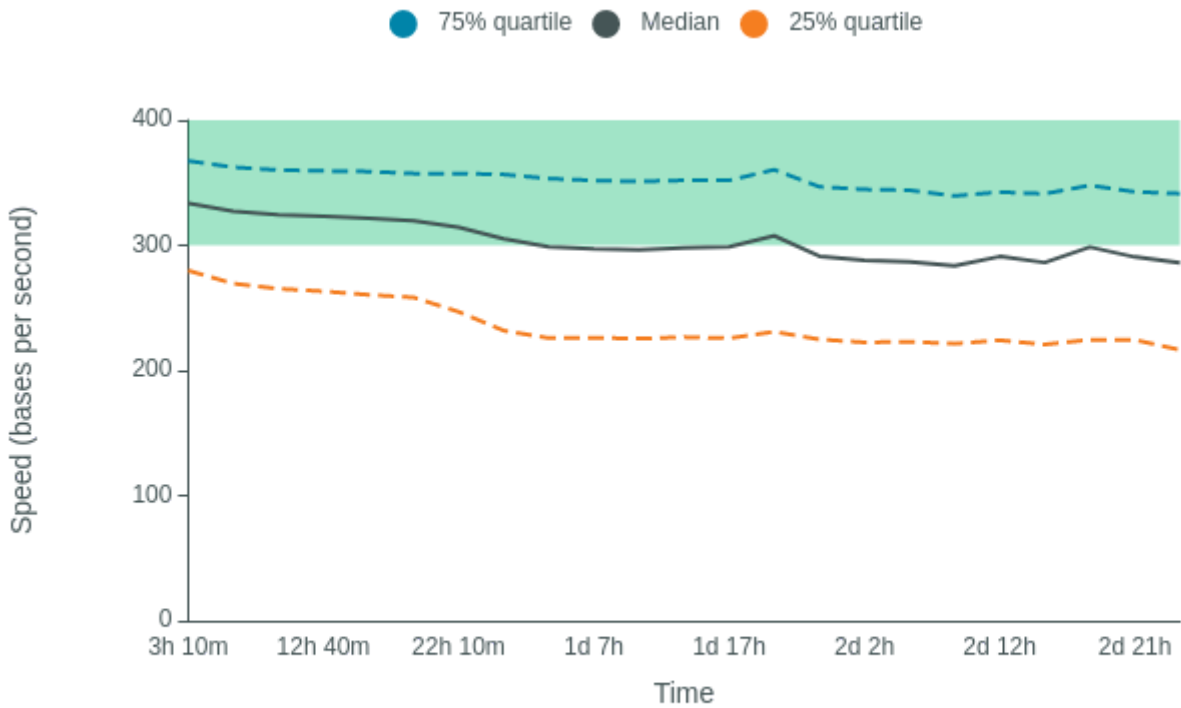


QScore

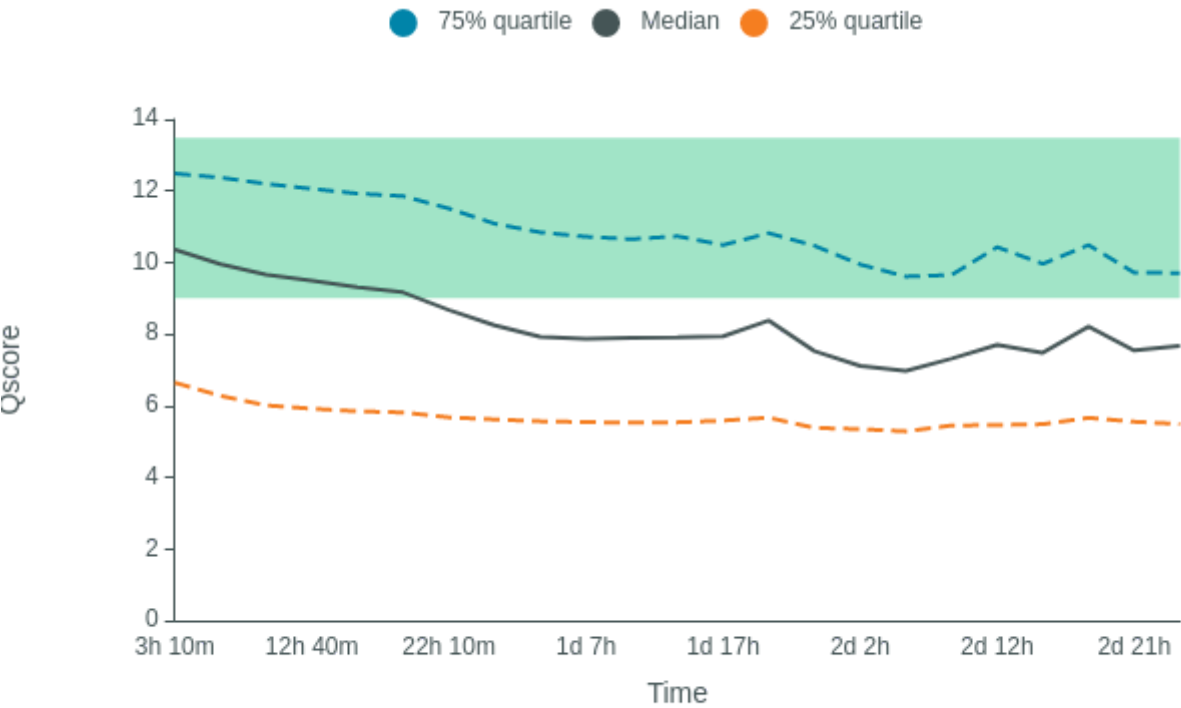



Disk Write Performance

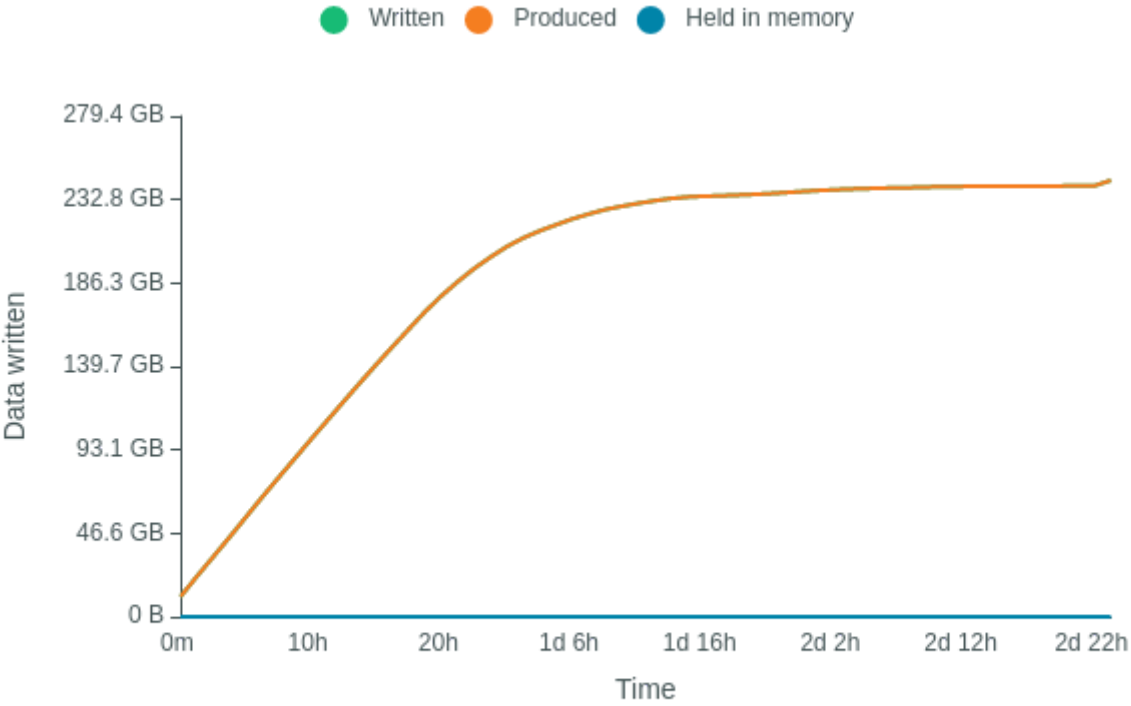

## **Run Debug Messages**

- The sequencing run has finished, but basecalling may continue July 25, 10:57
- Mux scan for flow cell FAL83665 has found a total of 24 pores. 24 pores available for immediate sequencing July 25, 09:57
- Performing Mux Scan July 25, 09:55
- Mux scan for flow cell FAL83665 has found a total of 28 pores. 26 pores available for immediate sequencing July 25, 08:25
- Performing Mux Scan July 25, 08:22
- Mux scan for flow cell FAL83665 has found a total of 36 pores. 35 pores available for immediate sequencing July 25, 06:52
- Performing Mux Scan July 25, 06:50
- Mux scan for flow cell FAL83665 has found a total of 29 pores. 28 pores available for immediate sequencing July 25, 05:20
- Performing Mux Scan July 25, 05:17
- Mux scan for flow cell FAL83665 has found a total of 29 pores. 28 pores available for immediate sequencing July 25, 03:47
- Performing Mux Scan July 25, 03:45
- Mux scan for flow cell FAL83665 has found a total of 36 pores. 33 pores available for immediate sequencing July 25, 02:15
- Performing Mux Scan July 25, 02:12
- Mux scan for flow cell FAL83665 has found a total of 50 pores. 50 pores available for immediate sequencing July 25, 00:42
- Performing Mux Scan July 25, 00:40
- Mux scan for flow cell FAL83665 has found a total of 61 pores. 58 pores available for immediate sequencing July 24, 23:10
- Performing Mux Scan July 24, 23:07
- Mux scan for flow cell FAL83665 has found a total of 70 pores. 67 pores available for immediate sequencing July 24, 21:37
- Performing Mux Scan July 24, 21:35
- Mux scan for flow cell FAL83665 has found a total of 62 pores. 53 pores available for immediate sequencing July 24, 20:05
- Performing Mux Scan July 24, 20:02
- Mux scan for flow cell FAL83665 has found a total of 63 pores. 61 pores available for immediate sequencing July 24, 18:32
- Performing Mux Scan July 24, 18:30
- Mux scan for flow cell FAL83665 has found a total of 64 pores. 60 pores available for immediate sequencing July 24, 17:00
- Performing Mux Scan July 24, 16:57
- Mux scan for flow cell FAL83665 has found a total of 65 pores. 60 pores available for immediate sequencing July 24, 15:27
- Performing Mux Scan July 24, 15:25
- Mux scan for flow cell FAL83665 has found a total of 75 pores. 68 pores available for immediate sequencing July 24, 13:55
- Performing Mux Scan July 24, 13:52
- Mux scan for flow cell FAL83665 has found a total of 70 pores. 63 pores available for immediate sequencing July 24, 12:22
- Performing Mux Scan July 24, 12:20
- Mux scan for flow cell FAL83665 has found a total of 73 pores. 62 pores available for immediate sequencing July 24, 10:50
- Performing Mux Scan July 24, 10:47
- Mux scan for flow cell FAL83665 has found a total of 73 pores. 63 pores available for immediate

- sequencing July 24, 09:17
- Performing Mux Scan July 24, 09:15
- Mux scan for flow cell FAL83665 has found a total of 71 pores. 65 pores available for immediate sequencing July 24, 07:45
- Performing Mux Scan July 24, 07:42
- Mux scan for flow cell FAL83665 has found a total of 67 pores. 61 pores available for immediate sequencing July 24, 06:12
- Performing Mux Scan July 24, 06:10
- Mux scan for flow cell FAL83665 has found a total of 79 pores. 71 pores available for immediate sequencing July 24, 04:39
- Performing Mux Scan July 24, 04:37
- Mux scan for flow cell FAL83665 has found a total of 107 pores. 96 pores available for immediate sequencing July 24, 03:07
- Performing Mux Scan July 24, 03:05
- Mux scan for flow cell FAL83665 has found a total of 163 pores. 139 pores available for immediate sequencing July 24, 01:34
- Performing Mux Scan July 24, 01:32
- Mux scan for flow cell FAL83665 has found a total of 180 pores. 157 pores available for immediate sequencing July 24, 00:02
- Performing Mux Scan July 23, 23:59
- Mux scan for flow cell FAL83665 has found a total of 186 pores. 154 pores available for immediate sequencing July 23, 22:29
- Performing Mux Scan July 23, 22:27
- Mux scan for flow cell FAL83665 has found a total of 218 pores. 173 pores available for immediate sequencing July 23, 20:57
- Performing Mux Scan July 23, 20:54
- Mux scan for flow cell FAL83665 has found a total of 277 pores. 207 pores available for immediate sequencing July 23, 19:24
- Performing Mux Scan July 23, 19:22
- Mux scan for flow cell FAL83665 has found a total of 346 pores. 235 pores available for immediate sequencing July 23, 17:52
- Performing Mux Scan July 23, 17:49
- Mux scan for flow cell FAL83665 has found a total of 357 pores. 230 pores available for immediate sequencing July 23, 16:19
- Performing Mux Scan July 23, 16:17
- Mux scan for flow cell FAL83665 has found a total of 412 pores. 262 pores available for immediate sequencing July 23, 14:47
- Performing Mux Scan July 23, 14:44
- Mux scan for flow cell FAL83665 has found a total of 537 pores. 319 pores available for immediate sequencing July 23, 13:14
- Performing Mux Scan July 23, 13:12
- Mux scan for flow cell FAL83665 has found a total of 678 pores. 373 pores available for immediate sequencing July 23, 11:41
- Performing Mux Scan July 23, 11:39
- Mux scan for flow cell FAL83665 has found a total of 747 pores. 359 pores available for immediate sequencing July 23, 10:09
- Performing Mux Scan July 23, 10:06
- Mux scan for flow cell FAL83665 has found a total of 823 pores. 380 pores available for immediate sequencing July 23, 08:36
- Performing Mux Scan July 23, 08:34
- Mux scan for flow cell FAL83665 has found a total of 926 pores. 410 pores available for immediate sequencing July 23, 07:04

- Performing Mux Scan July 23, 07:01
- Mux scan for flow cell FAL83665 has found a total of 960 pores. 407 pores available for immediate sequencing July 23, 05:31
- Performing Mux Scan July 23, 05:29
- Mux scan for flow cell FAL83665 has found a total of 980 pores. 410 pores available for immediate sequencing July 23, 03:58
- Performing Mux Scan July 23, 03:56
- Mux scan for flow cell FAL83665 has found a total of 1019 pores. 424 pores available for immediate sequencing July 23, 02:26
- Performing Mux Scan July 23, 02:23
- Mux scan for flow cell FAL83665 has found a total of 1033 pores. 417 pores available for immediate sequencing July 23, 00:53
- Performing Mux Scan July 23, 00:51
- Mux scan for flow cell FAL83665 has found a total of 1056 pores. 418 pores available for immediate sequencing July 22, 23:21
- Performing Mux Scan July 22, 23:18
- Mux scan for flow cell FAL83665 has found a total of 1084 pores. 427 pores available for immediate sequencing July 22, 21:48
- Performing Mux Scan July 22, 21:46
- Mux scan for flow cell FAL83665 has found a total of 1089 pores. 430 pores available for immediate sequencing July 22, 20:15
- Performing Mux Scan July 22, 20:13
- Mux scan for flow cell FAL83665 has found a total of 1118 pores. 432 pores available for immediate sequencing July 22, 18:43
- Performing Mux Scan July 22, 18:40
- Mux scan for flow cell FAL83665 has found a total of 1147 pores. 431 pores available for immediate sequencing July 22, 17:10
- Performing Mux Scan July 22, 17:08
- Mux scan for flow cell FAL83665 has found a total of 1146 pores. 434 pores available for immediate sequencing July 22, 15:38
- Performing Mux Scan July 22, 15:35
- Mux scan for flow cell FAL83665 has found a total of 1169 pores. 428 pores available for immediate sequencing July 22, 14:05
- Performing Mux Scan July 22, 14:02
- Mux scan for flow cell FAL83665 has found a total of 1185 pores. 433 pores available for immediate sequencing July 22, 12:32
- Performing Mux Scan July 22, 12:30
- Mux scan for flow cell FAL83665 has found a total of 1193 pores. 435 pores available for immediate sequencing July 22, 11:00
- Performing Mux Scan July 22, 10:57
- Starting sequencing procedure July 22, 10:57
- Waiting up to 300 seconds for temperature to stabilise at 34.0°C July 22, 10:53